


PREMIUM EXPERIENCE THE EXCELLENCE

PERFECT DONNING - due to special inner lining
HIGH WEARING COMFORT - due to especially developed shape
SAFE INSTRUMENT HANDLING - thanks to excellent grip

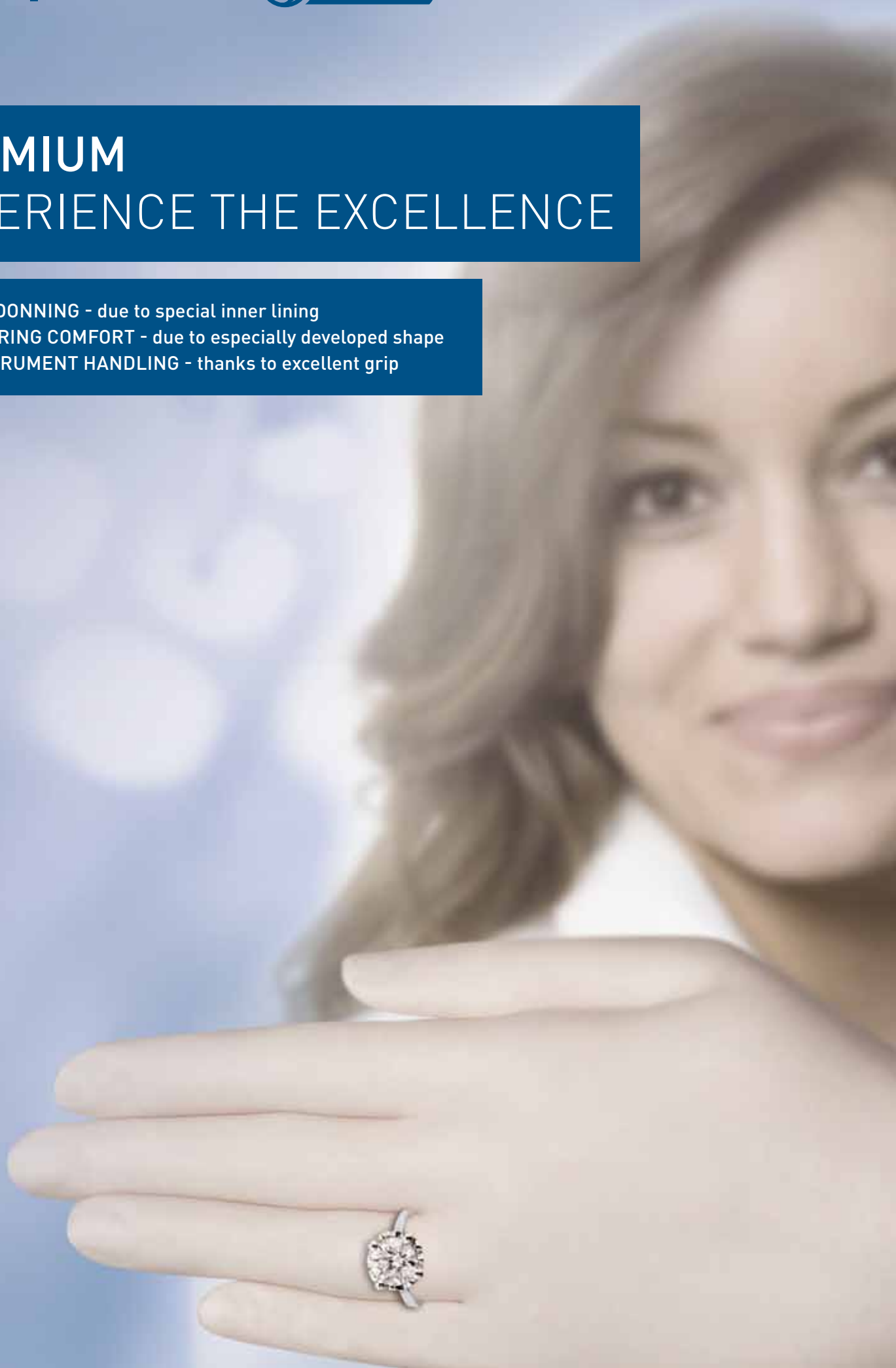


EXPERIENCE THE EXCELLENCE

For your professional top performance the best working material is just about good enough. Of course this credo is also true for our examination gloves. With the Sempercure® Premium we have developed a product which meets the highest expectations and sets new standards.

- For the Sempercure® Premium glove a new extra wide porcelain dipping mould with longer fingers was developed. The optimized length of the fingers significantly improves the **wearing comfort** of our powder free natural latex examination glove.
- Our powder free natural latex examination glove excels by its extraordinary, **well-gliding** surface. The textured fingers provide a good grip which is indispensable for **safe instrument handling**.
- The **synthetic inner lining** achieved by the special Sempermed technology makes it **easy** and **fast to don**.
- Like all products made by Sempermed, the Sempercure® Premium convinces with its particular **skin friendliness**.
- **The force at break¹⁾** is significantly higher than the minimum value required by the Standard and therefore affords the user a feeling of special protection.

¹⁾ The comparison refers to EN 455-2 (2011) and the average Sempercure® Premium production value, status 9/2011

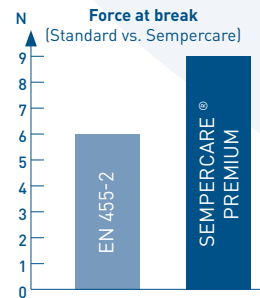

TECHNICAL DATA



sempercure® premium

TYPE	Unsterile, powder-free examination glove for single use
MATERIAL	Natural latex (NR)
COLOUR	Natural white
GLOVE SHAPE	Flat with rolled rim
OVERALL LENGTH/SIZE as per EN 455-2	XS, S, M, L, XL 240 mm
WALL THICKNESS imeasured in double layer in the palm area	median 0.26 mm
IMPERMEABILITY as per EN 455-1	AQL 1.5
FORCE AT BREAK as per EN 455-2 FORCE AT BREAK production value	≥ 6N median 9N
DURABILITY in original packaging if stored as per DIN 7716, ISO 2230	3 years
PACKAGING	Unsterile in dispenser carton: XS, S, M, L 100 pcs./XL 90 pcs. in transport carton: XS, S, M, L 1.000 pcs./XL 900 pcs./PU
MARKING	EN 1041, EN 980, EN 420 and ISO 15223
PURPOSE as set out in MDD 93/42 EEC as set out in PPE 89/686/EEC	Medical device class I Single-use protective glove - PPE category III (protection against chemical substances for limited time) ²⁾

CONTACT US!

Semperit Technische Produkte GmbH · Division Sempermed
Austria · A-1031 Vienna · Modecenterstraße 22 · Tel.: +43 1 79 777-0 · Fax: +43 1 79 777-630
e-mail: sempermed@semperit.at · www.sempermed.com